




DOUBLE WHEEL CUTTERS

Code N13xx

Sample cutter for strips of paper and board

Normas aplicables : ISO 1924 – SCAN P38

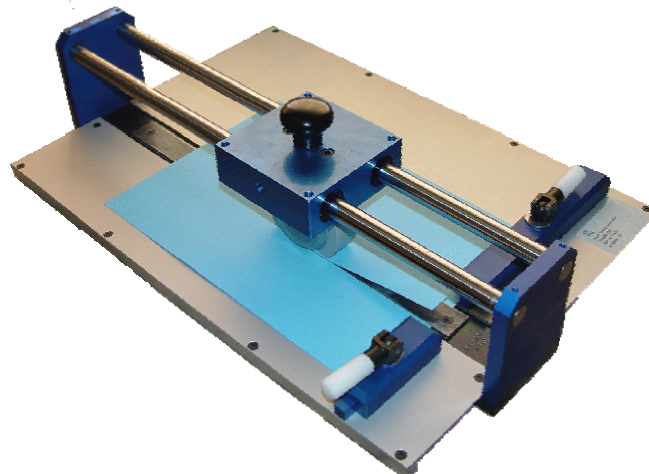

PTA SPAIN

C/ Gurutzegi Nº 16 Pol. Belartza
20018 SAN SEBASTIÁN
(Guipúzcoa) - Spain

PTA FRANCE

14, rue Evariste Galois / ZAC les Ruieres
38320 EYBENS - France

PTA ASIA

50 East Coast Road #02-73 Roxy Square
428769 SINGAPORE - Singapore



CHARACTERISTICS :

Strip cutter for the preparation of paper samples for tensile tests according to ISO 1924 and other standards.

The support with circular knives slides on two guide bars along a lower knife and cuts one (or optionally two) strips from the positioned paper sheet.

By default, the instrument delivered without fixing clamp. It can be ordered separately.

The instrument is available in different widths and cutting lengths.

Available cutting lengths : 350 or 550 mm

Available strip widths :

Ref 130½	width 12,7 mm
Ref 1315	width 15 mm
Ref 1325	width 25 mm
Ref 1301	width 25,4 mm
Ref 1330	width 30 mm
Ref 1350	width 50 mm
Ref 13100	width 100 mm
Ref 13135	width 135 mm
Ref 13150	width 150 mm
Ref 13xx	width > 8mm

WEIGHT AND DIMENSIONS (ONLY EQUIPMENT) :

290 x 470 x 270 mm (WxDxH) / 8 Kg

WEIGHT AND DIMENSIONS WITH PACKING BOX :

350 x 550 x 400 mm (WxDxH) / 15 Kg

PTA SPAIN

C/ Gurutzegi Nº 16 Pol. Belartza
20018 SAN SEBASTIÁN
(Guipúzcoa) - Spain

PTA FRANCE

14, rue Evariste Galois / ZAC les Ruires
38320 EYBENS - France

PTA ASIA

50 East Coast Road #02-73 Roxy Square
428769 SINGAPORE - Singapore